



Reprocessing Instructions for Reusable Surgical Instruments

English

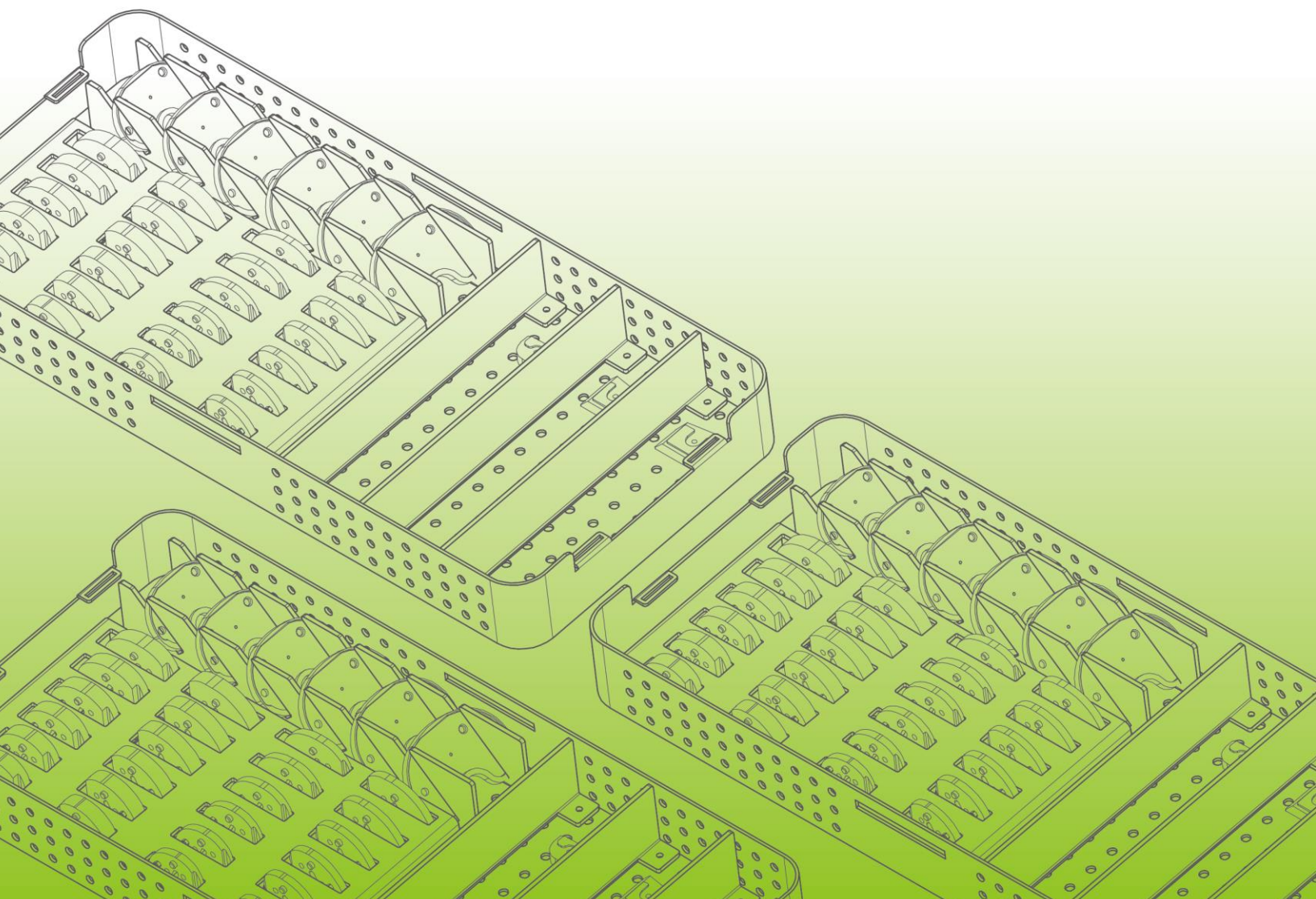



Table of contents

Description.....	1
Materials.....	1
 Cautions.....	2
Limitations on Reprocessing	3
Cleaning, Inspection, Packaging	4
A. Pre-treatment at the Point of Use	4
B. Containment and Transportation	4
C. Disassembly	4
D. Preparation of Cleaning (For All Cleaning Methods)	4
E. Cleaning- Manual Method with Ultrasonic	5
F. Cleaning -Automated /Mechanical Method	6
G. Inspection, Functional Check.....	6
H. Packaging	8
Sterilization	10
In the United States	10
Outside United States	10
Storage.....	11
Validation Information	11
For Manual Cleaning/Automated Cleaning Validation:	11
For Moist Heat Sterilization Validation:	11
Information.....	12
References	13
Symbols Glossary	14

DESCRIPTION

Reprocessing instructions are intended to assist hospital personnel in safe handling, and effective reprocessing **reusable surgical instruments**, instrument trays and cases supplied **Non-Sterile** by United Orthopedic Corporation (UOC). All hip, knee, trauma, and extremity of UOC reusable medical devices must be cleaning, inspection, packing and sterilization prior to use.

UOC has validated the processes provided in reprocessing instructions to be capable of being effective. Equipment, operators, cleaning agents and procedures all contribute to the efficacy of the processing. The healthcare facility should ensure that the selected reprocessing steps are safe and effective. Alternative methods of reprocessing outside the scope of these instructions may be suitable for reprocessing; however, those must be validated by the end user. In the event of conflicting national cleaning and sterilization requirements, these shall prevail over UOC recommendations.

MATERIALS

Aluminum alloy

Cobalt/Chromium alloys

Stainless steels

Titanium/Titanium alloys

Titanium Nitride coatings


Polymers

Silicone rubber



CAUTIONS

1. **Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner.**

Following symbol: 

2. Reprocessing instructions apply to:
 - **Reusable surgical instruments** supplied **Non-Sterile** by UOC
 - Instruments intended for reprocessing in a health care facility setting
3. **Do not** apply to single-use devices.
4. **Personal Protective Equipment (PPE)** should be worn when handling or working with contaminated or potentially contaminated devices. PPE includes: gown, mask, goggles or face shield and shoe covers.
5. **Do not** place heavy instruments on top of delicate devices.
6. **Do not** use metal brush, scouring pads, abrasive cleaner during manual cleaning procedures. These materials will damage the surface and finish of instruments. Soft-bristled, nylon brushes and pipe cleaner brushes should be used.
7. Cleaning agents with a pH 7~9, low foaming surfactants, nonabrasive, free-rinsing, biodegradable and environmentally friendly, provides for rapid soil dispersion or suspension should be used during cleaning to ensure that instruments are visible in the cleaning solution. Cleaning agents must be easily and completely rinsed from device surfaces to prevent accumulation of detergent residue.
8. Manual scrubbing with brushes should always be performed with the instrument below the surface of the cleaning solution to prevent formation of aerosols and splashing which may spread contaminants.
9. **Do not** use saline and cleaning/disinfection agents containing aldehyde, mercury, active chlorine, bromine, or iodide. These are corrosive and should not be used.
10. **Do not** place or soak instruments in Ringers Solution.
11. **Do not** use oil-lubricants. Because these may: 1) coat microorganisms; 2) prevent direct contact of the surface with steam; and 3) be difficult to remove.
12. **Do not** apply that the devices are not manufactured and/or distributed by UOC. Only devices manufactured and/or distributed by UOC should be placed in UOC instrument trays and cases. These validated reprocessing instructions are not applicable that the

- devices are not manufactured and/or distributed by UOC.
13. **Do not** use descaling agents that include morpholine in steam sterilizers. These agents leave residue which may damage polymer instruments over time.
 14. **Do not** stack cases or trays during sterilization. Because it may limit steam penetration and prevent effective sterilization of the instruments.
 15. Under certain classifications of risk, the World Health Organization (WHO), or local regulatory authorities recommend special Creutzfeldt-Jakob disease (CJD/TSE) inactivation processing procedures. Consult WHO and local regulations for further information.
 16. **Do not** be equal to/or greater than **140°C/284°F** in washer/sterilizers. The most polymers will occur severe surface damage.
 17. Flash (immediate-use) steam sterilization shall be used as an emergency procedure. Instruments shall be cleaned and disassembled without sterilization wrap or rigid container.
 18. May have sharp edges. Handle with care.
 19. For EU users, if any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

LIMITATIONS ON REPROCESSING

The useful life of these devices depends on many factors including the method and duration of each use, and the handling between uses. Careful inspection and functional test of the device before use is the best method of determining the end of serviceable life of the medical device. Evidence of damage and wear on a device may include but is not limited to corrosion (i.e. rust, pitting), discoloration, excessive scratches, flaking, wear and cracks. Improperly functioning devices, devices with unrecognizable markings, missing or removed (buffed off) part numbers, damaged, excessively worn and not repaired or reconditioned by UOC authorized personnel devices should not be used. If there is a need to return instruments to your sales or UOC, the devices must be cleaned, packaged and sterilized before returning.

CLEANING, INSPECTION, PACKAGING

<p>A. Pre-treatment at the Point of Use</p>	<ul style="list-style-type: none"> ● Remove contamination: As soon as possible after use, remove excessive soiling with a disposable wipe, rinse and flush lumened devices with sterile or deionized water to prevent the drying of soil and/or debris to the inside. ● Ensure that no instruments or pieces of instruments are left in the surgical site prior to closure, as they may not be detectable using imaging techniques such as X-ray or MRI and patient injury may result.
<p>B. Containment and Transportation</p>	<ul style="list-style-type: none"> ● Process instruments as soon as is reasonably possible after use. It is recommended not to delay cleaning for more than 2 hours. ● If transfer to the reprocessing area likely to be delayed, consider covering the medical devices with a damp cloth or store the medical devices in closed boxes to avoid drying of soil. ● If desired, place the instruments in its respective position within the instrument tray. The position of the instrument is labeled in its intended position within the tray.
<p>C. Disassembly</p>	<ul style="list-style-type: none"> ● Disassembly multi-piece or complex instruments referred to their cleaning instructions. Care should be executed to avoid losing small screws and components. The cleaning instructions are available from your local sales representative or UOC.
<p>D. Preparation of Cleaning (For All cleaning methods)</p>	<ul style="list-style-type: none"> ● Prepare a cleaning solution (Enzymatic/Detergent) proven efficacy and neutral pH 7~9 in accordance to the manufacturer's instructions provided information concerning specific materials, temperature, water quality, time, cleaning method. ● Tap Water could be referred to the utility water of AAMI TIR34: Hardness <150ppm, Conductivity <500µS/cm, Chlorides <250ppm, Total organic carbon <1mg/mL. Used primarily for flushing, washing and rinsing. The cold tap water is at less 40°C. The warm tap water is at 30~ 44°C. ● Deionized(DI) water (ultra-filter (UF), reverse osmosis (RO), or equivalent) shall be referred to the critical water of AAMI

TIR34[10], Hardness <1ppm, Conductivity <10µS/cm, Chlorides <1ppm, Bacteria <10CFU/mL, Endotoxin <10CFU/mL, Total organic carbon <1mg/mL. Used primarily for final rinses and steam generation. The warm deionized water is at 30~ 44°C. The Hot deionized water is at more than 44°C.

- Equipment: various sized soft-bristled brushes, lint-free cloths syringes, pipettes and/or water jet, ultrasonic cleaner, cleaning bath or vessel large enough to allow complete immersion of the instruments.

**E. Cleaning-
Manual
Method with
Ultrasonic**

- 1) **Soak** soiled instruments and prevent air bubbles to ensure that all surfaces have contact with an enzyme solution for a minimum recommended time specified by the enzymatic cleaning solution manufacturer or 20 minutes, whichever is longer.
- 2) **Brush** the instruments with cleaned soft-bristled, nylon brush to clean and remove all traces of blood and debris. Particular attention must be given to crevices, lumens, mated surfaces, connectors and other hard-to-clean areas. Lumens should be cleaned with a long, narrow, soft-bristled brush (i.e. pipe cleaner brush). For flexible shafts and springs, flex and relax the instrument under the cleaning solution while brushing.
- 3) **Flush** each difficult brush area thoroughly and aggressively in cold tap water for a minimum of 30 seconds. Use of a syringe or water jet will improve flushing of difficult to reach areas and closely mated surfaces. Repeat step 2~3 until all visible soil has been removed.
- 4) **Rinse** the instruments in cold top water for a minimum of 3 minutes.
- 5) **Ultrasonic cleaning:** Fully soak the opened instruments and prevent air bubbles to ensure that all surfaces have contact with an enzyme cleaning solution for 30 minutes at 45-50 kHz.
- 6) **Rinse** the instruments in cold tap water for a minimum of 3 minute. Repeat steps 2–6, until no visible debris, soil, enzyme cleaning solution remains on device.
- 7) **Final rinse** the instruments in **warm** deionized water for a minimum of 3 minutes to irrigate the challenging design

	<p>features.</p> <p>8) Dry the instruments after final rinse with a clean towel or clean compressed air until visibly dry.</p>																								
<p>F. Cleaning -Automated /Mechanical Method</p>	<p>1) Cleaning-Manual method with Ultrasonic step 1-6 should occur prior to this step. Mechanical washer with approved efficacy (e.g. Verified by ISO 15883) should be used.</p> <p>2) Load the opened instruments in the washer such can drain (for example, hinges should be open and cannulations and holes can drain).</p> <p>3) Process parameters follow "INSTRUMENT" cycle parameters validated by mechanical washer manufacturer and a pH neutral cleaning agent intended for use in automated cleaning using the MINIMUM cycle parameter set points below:</p> <table border="1" data-bbox="527 850 1448 1537"> <thead> <tr> <th>Cycle</th> <th>Minimum Time</th> <th>Temperature</th> <th>Liquid</th> </tr> </thead> <tbody> <tr> <td>Pre-wash</td> <td>2 minutes</td> <td><40°C</td> <td>Tap water</td> </tr> <tr> <td>Wash I</td> <td>5 minutes</td> <td>30~44°C</td> <td>Enzymatic solution</td> </tr> <tr> <td>Wash II</td> <td>5 minutes</td> <td>66°C</td> <td>Neutral pH Detergent solution</td> </tr> <tr> <td>Rinse</td> <td>10 minutes</td> <td>>44°C</td> <td>Deionized water</td> </tr> <tr> <td>Drying</td> <td>30 minutes</td> <td>100°C</td> <td>None</td> </tr> </tbody> </table>	Cycle	Minimum Time	Temperature	Liquid	Pre-wash	2 minutes	<40°C	Tap water	Wash I	5 minutes	30~44°C	Enzymatic solution	Wash II	5 minutes	66°C	Neutral pH Detergent solution	Rinse	10 minutes	>44°C	Deionized water	Drying	30 minutes	100°C	None
Cycle	Minimum Time	Temperature	Liquid																						
Pre-wash	2 minutes	<40°C	Tap water																						
Wash I	5 minutes	30~44°C	Enzymatic solution																						
Wash II	5 minutes	66°C	Neutral pH Detergent solution																						
Rinse	10 minutes	>44°C	Deionized water																						
Drying	30 minutes	100°C	None																						
<p>G. Inspection, Functional Check</p>	<ul style="list-style-type: none"> ● Visual inspect all instruments, instrument trays, cases after processing prior to sterilization for: (Generally un-magnified visual inspection under good light conditions.) <ul style="list-style-type: none"> ■ Cleanliness: Ensure the complete removal of soil from surfaces, tubes and holes, moveable parts. If soil is still present, re-clean the instrument. Particular attention should be paid to: Soil “traps” such as mating surfaces, 																								





hinges, shafts of flexible reamers; recessed features (holes, cannulations) ; features where soil may be pressed into contact with the device, e.g. drill flutes adjacent to the cutting tip, sides of teeth on broaches and rasps.


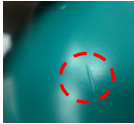






Note: If contamination is noted repeat the cleaning process.

- **Completeness, damage and/or excessive wear:** Visually inspect for no damage including but not limited to, malfunction, burrs, wear, tear, corrosion (rust, pitting), discoloration, creaked seals, excessive scratches and flaking.

Note: If damage or wear is noted that may compromise the function of the instrument, contact your sales or UOC.

- **Functional check** should be performed. Improperly functioning devices, devices with unrecognizable markings, missing or removed (buffed off) part numbers, damaged and worn devices should not be used. The disassembled instruments should be **reassembled to check the function.**

Check Items	Suitable for Use 	Not Suitable for Use 
Assemblable	Reassembled successfully.	Can't reassembled such as missing parts or can't smooth move.
Hinged instruments	Smooth operation throughout the intended range of motion.	Can't operate as intended.
Locking mechanisms	Can be locked to the position as intended and not loosening.	Unable to lock.
Cutting features Reamer/Drill	Edges should be continuous without burrs. 	Edge distortion/large nicks. 

Trials	Smooth and free of cracks and deep nicks. (Minor crack don't affect the smooth operation is acceptable.) 	Cracks and deep nicks. 
Hammering surfaces	Without burrs and large nicks.	Loose burrs or large nicks. 
Flexible features	Flex and relax to inspect flexible feature like shafts, springs for any damage and major deformation 	Major deformation or can't operate smoothly. 
Instrument tray/case	Inspect for burrs and locking from damage the wrap. (Minor burr is acceptable.) 	Burrs or locking mechanism loosening which may damage the wrap. 
UDI Information 	Use DataMatrix barcode scanners to scan the UDI carrier to check is readable.	Can't read UDI carrier.
<p>Note: If any status above is noted, do not use the instrument and contact your sales or UOC. The instruments can be reused only if it is reconditioned by UOC.</p>		

H. Packaging

- Use instrument cases to contain cleaned, checked, drying instruments that are provided in the instrument sets.
- For a sterilization wrap: UOC cases should be double wrapped according to AAMI/CSR technique. The packaging for terminally sterilized re-usable instruments should meet the following requirements: ISO 11607-1, CE Mark or FDA clearance, suitable for steam sterilization. **In the United States (US), only use an FDA-cleared Sterilization wrap.**
- For a rigid steam sterilization container system: Aesculap

SterilContainer System -JN-400~JN-446 included with base vents use as packaging for the instrument sets. No more than one case can be placed directly into a rigid steam sterilization container.

- **The total weight of a wrapped instrument case should not exceed 11.4kg/25lbs.**

STERILIZATION

Using ANSI/AAMI ST79 *Comprehensive Guide to Steam Sterilization and Sterility Assurance in Health Care* is recommended to be steam sterilized by the health care facility using the process parameters below, which has been validated by UOC under laboratory conditions to provide a 10⁻⁶ sterility assurance level (SAL). A verified, properly maintained and calibrated steam sterilizer is recommended. The process parameters of sterilization should be followed explicitly. It is the responsibility of the medical facility to ensure that reprocessing is performed using the appropriate equipment and materials, and that personnel in the reprocessing facility have been adequately trained.

- **In the United States**

Cycle Type	Pulse Number	Exposure Temperature	Minimum Exposure time	Minimum Drying times
Dynamic-Air-Removal (Pre-vacuum)	4	132°C (270°F)	4 minutes	30 minutes

- **Outside United States**

Cycle Type	Pulse Number	Exposure Temperature	Minimum Exposure time	Minimum Drying times
Dynamic-Air-Removal (Pre-vacuum)	4	132~135°C (270~275°F)	4 minutes	30 minutes
Dynamic-Air-Removal (Pre-vacuum) For CJD/TSE contamination ¹		134°C ¹ (273°F)	18 minutes ¹	

¹ Under certain classifications of risk, the World Health Organization (WHO), or local regulatory authorities recommend special Creutzfeldt-Jakob disease (CJD/TSE) inactivation processing procedures. Consult WHO and local regulations for further information.

STORAGE

After sterilization, re-usable instruments should be stored the in the sterilization wrap in a dry and dust-free place. The shelf life is depending on the sterile barrier employed, storage manner, environmental conditions and handling. A maximum shelf life for sterilized re-usable instruments before use should be defined by each health care facility.

VALIDATION INFORMATION

● For Manual Cleaning/Automated Cleaning Validation:

The reusable devices are exposed to the artificial test soil which is a mixture of Edinburgh test soil and bone meal to simulate the clinical soil, and be dried for 2 hours prior to clean. After manual cleaning/automated cleaning, visually inspect for any sign of remaining soil and test residual protein and total organic carbon (TOC) to an acceptable level.

- Following standards: AAMI TIR12, AAMI TIR30, AAMI TIR34, ANSI/AAMI ST79, ANSI/AAMI ST81, ISO 17664 and "Reprocessing Medical Devices in Health Care Settings: Validation Methods and Labeling Guidance for Food and Drug Administration Staff".
- Detergents: 3M Rapid Multi-Enzyme Cleaner(70503), Prolystica 2X Concentrate Enzymatic Presoak And Cleaner (1C33), Prolystica 2X Concentrate Neutral Detergent(1C32)
- Automated washer – Miele PG 8535

● For Moist Heat Sterilization Validation:

The reusable devices in a maximum weight of loaded tray are placed biological indicators (BI) inoculated with more than one million (10^6) resistant spores (Geobacillus stearothermophilus) in the most challenging locations. The tray is packed with a sterilization wrap or rigid steam sterilization container system. The packed tray loaded in an empty chamber of steam sterilizer is validated a 10^{-6} sterility assurance level (SAL) by overkill method that half-cycle sterilization cycle result in total kill of all BIs demonstrating in accordance with ISO 17665-1 Annex D. Recommended dry time is validated by demonstrating that pre-sterilization and post sterilization wrap/filter weight is no exceed ± 3 percent weight gain and no visible moisture on or within the sterilized packing, container, tray, instruments.

Following standards: AAMI TIR12, ANSI/AAMI ST77, ANSI/AAMI ST79, ANSI/AAMI ST81, ISO 17664, BS EN ISO 17665-1, ISO/TS 17665-2, AAMI/ISO TIR17665-3 and "Reprocessing Medical Devices in Health Care Settings: Validation Methods and Labeling Guidance for Food and Drug Administration Staff".

- Sterilization wrap: Kinguard KC600 One-Step wrap
- Rigid steam sterilization container system: Aesculap SterilContainer Systems JN446(Bottom), JK486(Lid), MD344(Filter).
- Steam Sterilizer: AMSCO Eagle 3023-S Vacumatic Prevacuum Steam Sterilizer.

INFORMATION

For further information, please contact



United Orthopedic Corporation
No. 57, Park Ave 2, Science Park,
Hsinchu City, 30075, Taiwan.
TEL: +886-3-5773351
FAX: +886-3-5777156
E-MAIL: service@uoc.com.tw












mdi Europa GmbH
Langenhagener Straße 71
30855 Langenhagen, Germany
TEL: +49-511-3908 9530
FAX: +49-511-3908 9539

REFERENCES

1. AAMI TIR12: Designing, testing, and labeling reusable medical devices for re-processing in health care facilities: A guide for medical device manufacturers.
2. AAMI TIR30: A compendium of processes, materials, test methods, and acceptance criteria for cleaning reusable medical devices.
3. AAMI TIR34: Water for the reprocessing of medical devices
4. ANSI/AAMI ST77: Containment devices for reusable medical device sterilization.
5. ANSI/AAMI ST79: Comprehensive guide to steam sterilization and sterility assurance in health care facilities.
6. ANSI/AAMI ST81: Sterilization of medical devices—Information to be provided by the manufacturer for the processing of re-sterilizable medical devices.
7. ASTM F565: Standard Practice for Care and Handling of Orthopedic Implants and Instruments.
8. ISO 11607-1: Packaging for terminally sterilized re-usable instruments.
9. ISO 15223-1: Medical devices- Symbols to be used with medical device labels, labeling and information to be supplied - Part 1: General requirements.
10. ISO 17664: Sterilization of re-usable instruments -Information to be provided by the manufacturer for the processing of re-sterilizable re-usable instruments.
11. ISO 17665-1: Sterilization of health care products - Moist heat - Part 1: Requirements for the development, validation and routine control of a sterilization process for medical devices.
12. Reprocessing Medical Devices in Health Care Settings: Validation Methods and Labeling Guidance for Food and Drug Administration Staff, 2015.

SYMBOLS GLOSSARY

Symbol(s)	Title of symbol	Description of symbol	Standard / Ref. No.
	Medical Device	Instrument, apparatus, implement, machine, appliance, implant, reagent for <i>in vitro</i> use, software, material or other similar or related article, intended by the manufacturer to be used, alone or in combination, for human beings, for one or more of the specific medical purpose.	ISO 15223-1
	Catalogue number	Indicate the manufacturer's catalogue number.	ISO 15223-1 / 5.1.6
	Batch code	Indicate the manufacturer's batch code.	ISO 15223-1 / 5.1.5
	Manufacturer	Indicate the medical device manufacturer, as defined in EU Directives 93/42/EEC.	ISO 15223-1 / 5.1.1
	Caution	Indicate the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.	ISO 15223-1 / 5.4.4
	Authorized representative in the European Community	Indicate the authorized representative in the European Community.	ISO 15223-1 / 5.1.2
	Rx only	Caution: Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner.	
Non-Sterile	Non-Sterile	Indicate the device is supplied non-sterile.	
	CE Marking	Indicate a presumption of conformity with the relevant EU REGULATION.	Directive 93/42/EEC
	CE Marking	Indicate a presumption of conformity with the relevant EU Directives. This symbol is applicable for Class IIa devices.	Directive 93/42/EEC
QTY	Quantity	Indicate the quantity in a package of the medical device.	



Each Step **We Care**

Contact Us

■ Taiwan

United Orthopedic Corporation

12F, No.80, Sec.1, Chenggong Rd., Yonghe Dist., New Taipei City 23452, Taiwan

Tel: +886 2 2929 4567 Fax: +886 2 2922 4567

■ China

United Medical Instrument Co., Ltd.

Room 402, Building 7, No.697, Ling Shi Road, ShangHai 200072, China

Tel: +86 21 3661 8055 Fax: +86 21 6628 3031

■ Japan

United Orthopedic Japan Inc.

Ginyo Bldg 4th floor, 2-9-40 Kitasaiwai, Nishi-Ku Yokohama City, Kanagawa, 220-0004, Japan

Tel: +81 45 620 0741 Fax: +81 45 620 0742

■ USA

UOC USA Inc.

20, Fairbanks, Suite 173, Irvine 92618, CA
Tel: +1 949 328 3366 Fax: +1 949 328 3367

■ EU Representative

mdi Europa GmbH

Langenhagener Strasse 71, Langenhagen 30855, Germany

Tel: +49 511 3908 9530 Fax: +49 511 3908 9539

■ Europe

UOC (Suisse) SA

Avenue Général Guisan 60A, Pully 1009, Switzerland

Tel: +41 21 634 70 70 Fax: +41 21 634 70 76

UOC (France)

21, Rue de la Ravinelle, Nancy 54000, France

Tel: +33 3 83 23 39 72 Fax: +33 3 83 23 39 10



www.uoc.com.tw

CE CE 2797

Copyright © 2020 United Orthopedic Corp. All rights reserved.